


# Introducing VektorPRO – Class II correction simplified.

The **ONLY** appliance that utilizes gentle, intrusive force vectors –  
just what is needed for efficient Class II correction.

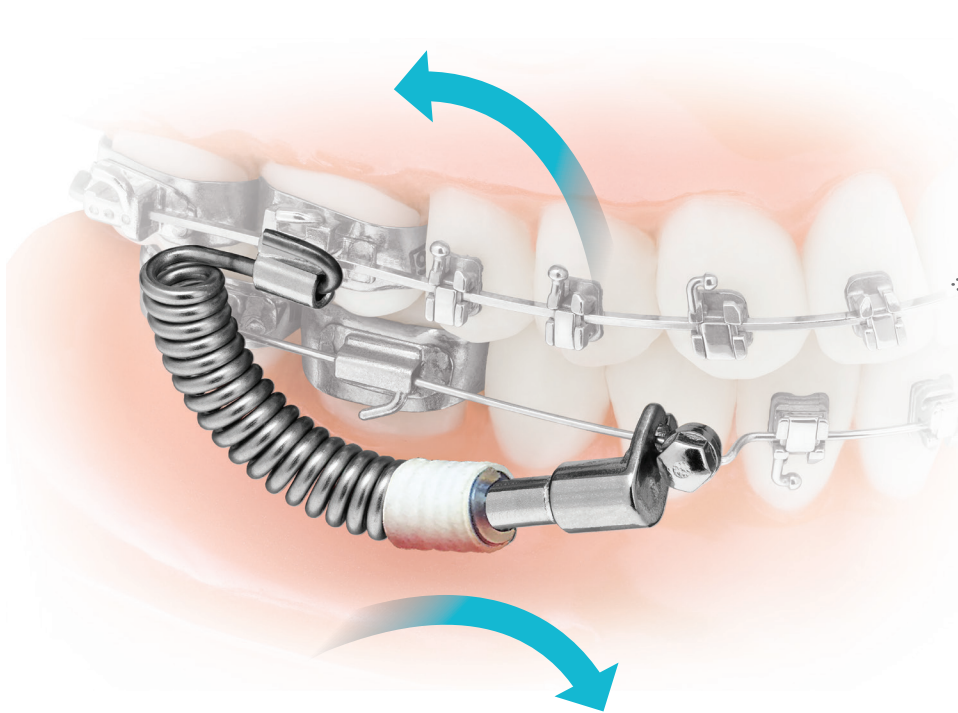

- The appliance flexes in **25-45% of its overall length, producing gentle, curved force vectors**
- Protective silicone sleeve prevents soft tissue irritation
- Robust nickel titanium to stainless steel assembly
- Easy installation
- Free from the “food zone”
- Allows natural lateral mandibular movement for patient comfort

Other appliances produce extrusive tipping forces on the maxilla and anterior teeth, which is always contraindicated for the over erupted Class II patient. Since the vector control module (VCM) is made from nickel titanium, the appliance is break-resistant and force values of only 3.5 ounces are produced. Provide your patients the latest Class II correction alternative with the VektorPRO appliance.


See why the VektorPRO appliance is in a  
class of its own. Call or visit us online:

800-348-8856 | [www.tportho.com](http://www.tportho.com)

# VektorPRO Appliance

The VektorPRO Class II Correction Appliance is a fixed intraoral device that achieves class II correction easily and efficiently with low forces. The patented design of the appliance introduces gentle intrusive vector forces - unlike other class II correctors that place extrusive vector forces on the anterior teeth.

The appliance produces gentle sweeping forces, and the mandible is free to move laterally resulting in patient comfort. The appliance is easy to install, clears the “food zone” and is compliance-independent.

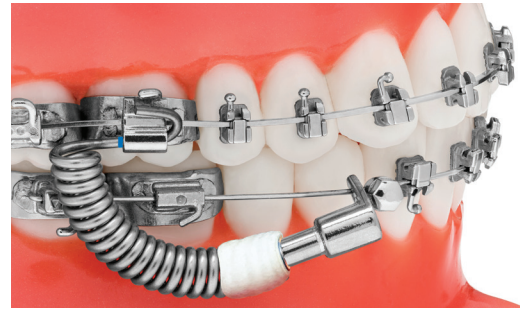

Patented

<b>Kits</b>	<b>Order No.</b>
Class II Corrector Multi-Kit Kits consists of 2 pairs size 27, 2 pairs size 30, 2 pairs size 33, 1 archwire installation tool, 20 archwire locks	<b>610-500P</b>
Class II Corrector Kit—size 27	<b>610-527P</b>
Class II Corrector Kit—size 30	<b>610-530P</b>
Class II Corrector Kit—size 33	<b>610-533P</b>

Single Kit includes Left / Right Component plus 4 Archwire Locks



Multi-Kit

<b>Component Part</b>	<b>Order No.</b>
27mm Left	<b>610-527PL</b>
27mm Right	<b>610-527PR</b>
30mm Left	<b>610-530PL</b>
30mm Right	<b>610-530PR</b>
33mm Left	<b>610-533PL</b>
33mm Right	<b>610-533PR</b>
Package of 5	
Archwire Locks	<b>610-502</b>
Package of 10	



Left / Right Component



Archwire Lock

<b>Tool</b>	<b>Order No.</b>
Archwire Lock Installation Tool	<b>610-501</b>
Each	



Archwire Lock  
Installation Tool

