

Featured Products from TPO®

INNOVATIONS IN ORTHODONTICS



TP Orthodontics, Inc.

Trusted
innovative
products
and
solutions.

Note: All TPO products may not be available in your area.
Contact your TPO representative for more information.



ORTHODONTIC INNOVATION Starts Here

TP Orthodontics offers some of the most innovative products and solutions available. From the original Tooth Positioner that inspired the concept of aligner treatment, to the latest in aesthetic brackets featuring *Personalized Color-Matching Technology* brand brackets to exceptional finishing solutions, we strive to develop unique solution-based concepts in orthodontics.

Trust That's Been Earned.

When you choose TP Orthodontics, you can be assured the products used in your practice are designed for efficiency, dependability and overall patient satisfaction. That's how we built our reputation as one of the most trusted suppliers of quality orthodontic products and solutions – and the reason we've been at this since 1942.

Unmatched Quality Standards.

TPO® maintains a worldwide Quality Management System that adheres to international standards and governmental regulations required for medical devices. This means you can be confident the products you receive from TP Orthodontics are held to the most rigorous standards in the industry – not every supplier can say that!

As the most experienced manufacturer in the industry, TP Orthodontics is a name you can trust.



AESTHETIC BRACKET SYSTEMS

The clear match for you AND your patients.

TPO provides you and your patients the ultimate combination of beauty and performance with our extensive line of aesthetic bracket systems.

Polychrystalline Ceramic Brackets

The doctor-preferred design coupled with patient preferred aesthetics.

- Patented, flexible polymer mesh base serves as a protective barrier between the ceramic bracket and enamel
- No special tools are required for debonding and no need to fracture the bracket
- Designed for multiple, secure ligation options



ClearVu® Mini Brackets
Roth or MBT



InVu® Mini Brackets
Roth or MBT



Tip-Edge® PLUS
Aesthetic Brackets

Virtually Invisible
Flexible Polymer Base
True-Twin Design
Pre-Coated
Easy Debonding
Multiple Ligation Options

Personalized Color-Matching Technology®

- Creates a look of transparency that matches virtually any tooth color
- Proprietary ceramic material closely duplicates optical properties of natural teeth
- Absorbs certain wavelengths of light, allowing only natural tooth color to show



Meet Patient Demand and Set Your Practice Apart

More patients than ever are requesting aligner treatment. *Personalized Color-Matching Technology* brand brackets from TPO are an excellent solution for patients who are not suitable candidates for aligner therapy. Use TPO ceramic brackets to start and keep your aesthetic cases on track – with predictable results. This unique solution is available with the following TPO bracket lines:

- ClearVu Cosmetic Brackets
- InVu Aesthetic Brackets
- Tip-Edge PLUS Aesthetic Brackets

Aesthetic Archwire

- Perfect for use with ceramic brackets
- Features a thin tooth-colored coating on the labial surface that does not alter the wire size
- Uncoated near the archwire slot to glide smoothly with no increase in friction
- Available in a variety of arch shapes and wire types
- Exclusive Convenience-Pak® hygienic packaging



METAL BRACKET SYSTEMS

Nu-Edge® Self-Ligating Brackets


The proven CoCr performance of Nu-Edge brackets with a self-ligating design.

- Cobalt chromium is biocompatible for nickel sensitive patients
- Provides audible assurance that the archwire is locked into place
- Maze base is designed for optimal bond strength
- Multiple ligation options

Nu-Edge Twin Brackets

Stronger CoCr alloy material. Lower friction. Superior compatibility.

- 99.9% nickel-free for sensitive patients
- Cobalt chromium design reduces friction forces
- Brackets are cast rather than machined or milled to enhance low-friction properties
- Dependable twin design

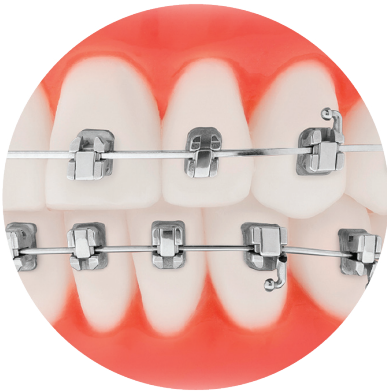

Nu-Edge Self-Ligating
Brackets



Nu-Edge Twin Brackets



Tip-Edge Metal Brackets

Tip-Edge® PLUS Brackets

Tip-Edge PLUS brackets are part of a versatile system that utilizes light forces to move teeth naturally and efficiently.


- Efficient treatment without the need for headgear, implants or functional appliances
- Pre-adjusted to capture 100% torque and tip – zero-tolerance finish
- Extremely light forces – only 1 to 3 ounces
- Automatic rotational control
- Eliminates unwanted torque reaction to adjacent teeth

Pre-Finisher® SEMI-CUSTOM FINISHING APPLIANCE

Off-the-shelf convenience, semi-custom finishing solution.

The Pre-Finisher appliance is a semi-custom solution perfect for improving the occlusion in aligner and fixed-appliance patients. When used as an adjunct to aligner treatment, it can reduce or eliminate additional trays required for case refinement – while correcting the posterior open bite often seen with aligner patients.

- No laboratory wait or fees; immediate delivery to the patient
- Relates teeth in both arches simultaneously; corrects open bite
- Aligns, rotates and closes spaces
- Sized in one-millimeter increments for precise fitting
- Optional seating springs facilitate placement and retention
- Ideal solution for the already-compliant aligner patient
- Results in 2-3 weeks

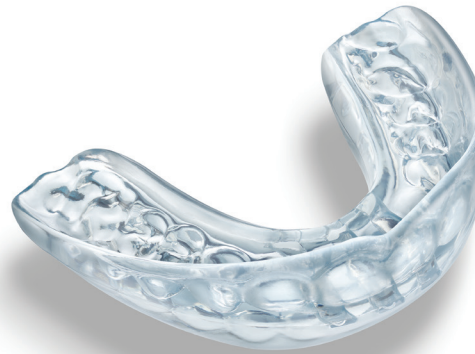

Inventory kit

Accessories

Crimpable Hooks



- The first slide-proof crimpable hook for rectangular wire to accept force loads equal to welded hooks
- Available in oriented packaging for ease of use
- **Coated and Ribbed Crimpable Hooks are also available with Aesthetic Coating.**



Archwires

- Wires available with or without aesthetic coating
 - Reflex® Nickel Titanium
 - Shiny Bright Stainless Steel
 - Braided 8-Strand Reflex
- Round and rectangular wire sizes
- Other wire options available



Ultra Low Profile Molar Tubes

- Extremely low profile for greater patient comfort
- Minimizes potential for occlusal force interference
- Funnel-shaped entrance for ease of wire insertion



TP Orthodontics, Inc.
100 Center Plaza
La Porte, IN 46350-9672 USA
800-348-8856

 tportho.com

 facebook.com/tportho

 twitter.com/tportho

 youtube.com/user/tportho

 linkedin.com/company/tp-orthodontics

Visit **tportho.com**
or contact one of our offices
below to learn more:

MAIN OFFICE

[Europe](#)
TP Orthodontics Europe
Fountain Court
12 Bruntcliffe Way
Morley, Leeds, LS27 0JG
England

Phone: 44 (0) 113 2539192
Fax: 44 (0) 113 2539193
E-mail: tpeng@tportho.com

EUROPEAN OFFICES

[Austria](#)
Phone/Fax: 0044 113 25 26 247

[France](#)
Phone/Fax: 0800 91 36 87

[Germany](#)
Phone/Fax: 0800 181 47 19

[Portugal](#)
Phone/Fax: 800 844 040

[Spain](#)
Phone/Fax: 930 17 32 57

[UK/Ireland](#)
Tel/Fax: 0044 1132539192

RESELLERS

A list of authorized resellers
can be found on our website at:
tportho.com/about/distributors



100-014ENG, Rev. 2 - 2/20

ClearVu, Convenience Pak, InVu, Nu-Edge, Personalized Color-Matching Technology, Pre-Finisher, Reflex, Tip-Edge and TPO are registered trademarks of TP Orthodontics, Inc.
ClearVu is manufactured under US Patents 6,685,468, 7,780,442; Brazil BR1120150234003; China ZL03824889.1; Europe 1352617, 1562511, 14775271.1; Germany 60315212.0, 60332263.8-08; Japan 4503440; Mexico 274623. InVu is manufactured under US Patents 6,685,468, 6,786,720, 7,780,442; Brazil BR1120150234003; Britain 1562511; China 03822189.6, ZL03824889.1; Europe 1352617, 14775271.1, 1545380; France 1562511; Germany 60332263.8-08; Italy 1562511; Japan 4414334, 4503440; Mexico 265725, 274623; Spain 1562511, ES2362089T3. Tip-Edge PLUS is manufactured under US Patents 6,682,345, 6,685,468, 6,786,720, 7,780,442; :Brazil BR1120150234003; China 03806949.0; Europe 14775271.1; Germany 60327861.2; Italy 1539017; Japan 4391242; Mexico 269871; Spain 1539017.

© 2020 TP Orthodontics, Inc. All rights reserved.