



TP Orthodontics, Inc.

[tportho.com](http://tportho.com)

[cosmeticbraces.com](http://cosmeticbraces.com)

# CLASS II CORRECTION


# CLASS II CORRECTION

**VektorPRO Class II Corrector**

**Flip-Lock® HERBST® Appliance**

# VektorPRO Class II Corrector

- Fixed intraoral device that achieves class II correction with low forces.
- Gentle intrusive vector forces – only 3.5 ounces are produced.
- Mandible is free to move laterally resulting in patient comfort.
- Easy to install, clears the “food zone” and is compliance-independent.

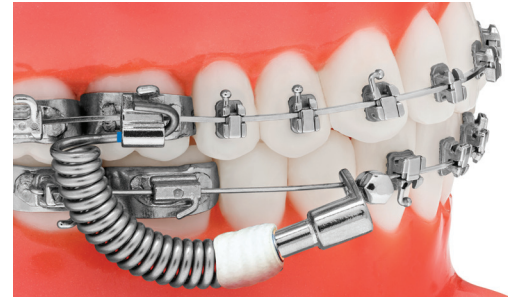

Patented

<b>Kits</b>	<b>Order No.</b>
Class II Corrector Multi-Kit Kits consists of 2 pairs size 27, 2 pairs size 30, 2 pairs size 33, 1 archwire installation tool, 20 archwire locks	<b>610-500P</b>
Class II Corrector Kit—size 27	<b>610-527P</b>
Class II Corrector Kit—size 30	<b>610-530P</b>
Class II Corrector Kit—size 33	<b>610-533P</b>
Single Kit includes Left / Right Component plus 4 Archwire Locks	



Multi-Kit

<b>Component Part</b>	<b>Order No.</b>
27mm Left	<b>610-527PL</b>
27mm Right	<b>610-527PR</b>
30mm Left	<b>610-530PL</b>
30mm Right	<b>610-530PR</b>
33mm Left	<b>610-533PL</b>
33mm Right	<b>610-533PR</b>
Package of 5	



Left / Right Component

<b>Archwire Locks</b>	<b>610-502</b>
Package of 10	



Archwire Lock

<b>Tool</b>	<b>Order No.</b>
Archwire Lock Installation Tool	<b>610-501</b>
Each	



Archwire Lock Installation Tool

CLASS II CORRECTION

# Flip-Lock® HERBST® Appliance

- Activates using the “Flip-Lock” mechanism.
- Place crimpable spacers at regular intervals onto the exposed piston.
- No pins or screws.
- Available with and without a cantilever arm.
- Optional archwire slots allow the appliance to be placed in conjunction with brackets.




Single Patient and Inventory Kits	Order No.
Single Patient Appliance Kit 2 pistons, 2 tubes, 4 ball hinges on small pad without slot	<b>610-400</b>
Single Patient Appliance Kit 2 pistons, 2 tubes, 4 ball hinges on large pad with slot	<b>610-410</b>
Single Patient Appliance Kit 2 pistons, 2 tubes, 4 ball hinges on large pad without slot	<b>610-430</b>
Single Patient Appliance Kit 2 tubes, 2 pistons, 2 posterior ball hinges, 2 cantilever/ball hinges without slot	<b>610-420</b>
Single Patient Appliance Kit 2 tubes, 2 pistons, 2 posterior ball hinges, 2 cantilever/ball hinges with slot	<b>610-424</b>
Inventory Kit 30 pistons, 30 tubes, 60 ball hinges, 95 spacers	<b>610-409</b>



Component Part	Order No.
Pistons	<b>610-403</b>
Tubes (5 left, 5 right)	<b>610-404</b>
Ball Hinge with small pad without slot	<b>610-405</b>
Posterior Ball Hinge without slot and large welding pad	<b>610-422</b>
Posterior Ball Hinge with slot and large welding pad	<b>610-423</b>
Cantilever/Anterior Ball Assembly without slot (5 left, 5 right)	<b>610-421</b>
One-piece Cantilever/Anterior Ball Assembly with slot (5 left, 5 right)	<b>610-425</b>

Package of 10

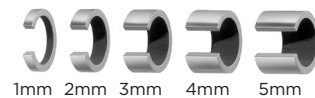

Crimpable Spacers	Order No.
1mm	<b>610-401</b>
2mm	<b>610-402</b>
3mm	<b>610-406</b>
4mm	<b>610-407</b>
5mm	<b>610-408</b>

Package of 10